

MICRONVIEW

BAS-Eco

BioAerosol Sampler-Economic



- Economic and portable. Designed for use in ISO and GMP standard cleanrooms and aseptic manufacturing areas
- 28.3, 100 and 200 LPM flow rate pre-selectable
- Complies with ISO 14698-1 and meets EU GMP Annex 1

Product Features



Accurate; Reliable; Economic

- Sealed anodized aluminum enclosure protects against moisture, dust and chemical ingress, and is resistant to disinfectants commonly used in cleanrooms.
- Air is impacted on the agar surface through the perforated cover; meeting isokinetic sampling requirements for a range of flow rates.
- In the event that there's a change in inlet pressure or resistance, a built-in high precision mass flowmeter will detect the air flow automatically and adjust the blower speed in real-time with millisecond frequency to ensure $\pm 2.5\%$ (28.3L/min, 100L/min) and $\pm 3\%$ (200L/min) flow accuracy.
- Reduced weight, compact screen design, meeting the clients' requirements for economic and efficient sampling.

Sampling

Distinct and Direct

- Quick set sampling specifications
- 4" Rapid response touch screen
- Supports remote sampling

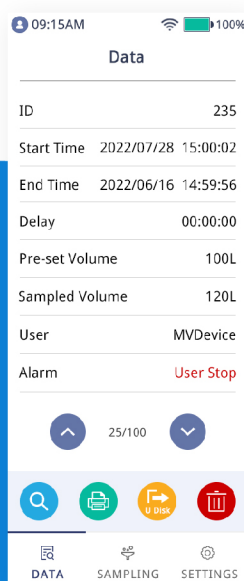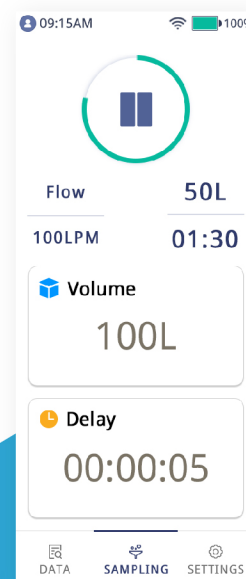

Data

Simple and Flexible


- Simple interface design, quick search capability
- One-click report export and printing
- Wireless connection with the APC and BAMS to support triggered sampling

Isokinetic Sampling

When performing isokinetic sampling, the sampler is designed so that the airflow direction into the inlet of the sampling head is consistent with the average velocity of the unidirectional flow through the sampling head.





The flow velocity above the sampling head is 0.4m/s, consistent with isokinetic sampling.

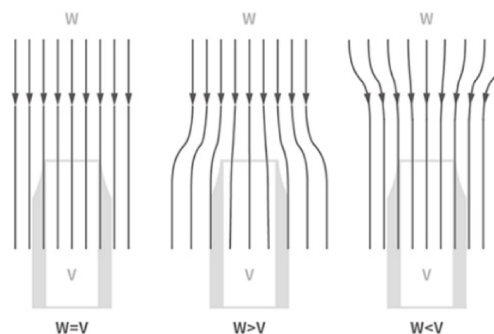

V=Sampling speed of sampling head
W=Speed of unidirectional flow

Remote Sampling

Used for sampling in an isolator, biological safety cabinet, horizontal laminar flow clean workbench and other significant zones.

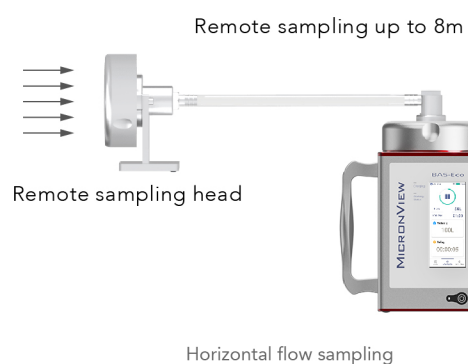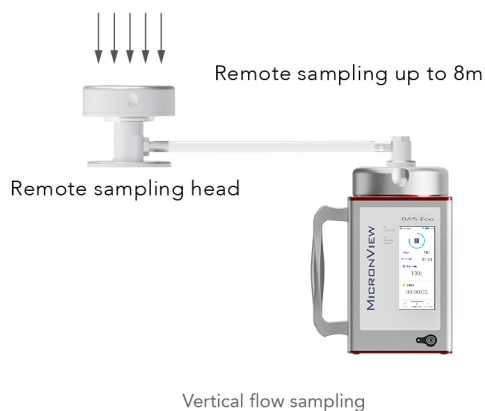

Compressed Gas Sampling

Compressed gas is decompressed to atmospheric pressure by the High Pressure Diffuser (HPD) and transported to the BioAerosol Sampler-Economic for sampling.



Specification Sheet

Specification	BioAerosol Sampler-Economic I BAS-Eco
Flow	28.3L/min $\pm$ 2.5%; 100L/min $\pm$ 2.5%; 200L/min $\pm$ 3%
Flow rate control	Built-in mass flow meter
Sampling volumes	1-6000L
Delay	0-24 hours
Sampling head	Φ 0.32mm*280 holes (28.3LPM) Φ 0.60mm*280 holes (100LPM) Φ 0.70mm*400 holes (200LPM) Anodized Aluminum(316L stainless steel available, customizable sampling head)
D50	0.79 μ m (28.3LPM; Φ 0.32mm*280 holes); 1.08 μ m (100LPM; Φ 0.60mm*280 holes); 1.15 μ m (200LPM; Φ 0.70mm*400 holes)
Weight / Dimensions	1.800Kg(3.968lbs)(Include Battery) Dimensions(without footmat) 218(H) $\times$ 120(W) $\times$ 120(D)mm; 8.58(H) $\times$ 4.72(W) $\times$ 4.72(D)inches
External surface	Anodized Aluminum
Power	AC 100-240V, 50 Hz/60 Hz
Battery	10.8V, 6800mAh rechargeable lithium battery Battery last $\geq$ 10 hours (28.3LPM) Battery last $\geq$ 10 hours (100LPM) Battery last $\geq$ 6 hours (200LPM)
Communication	USB, WIFI
Application extension	Wireless connection with the APC and BAMS to support triggerd sampling
Data storage	20,000 sets of records
Alarm	LED lighting and beep
Display	4" capacitive touch screen
Language	English, Chinese
Calibration frequency	Annual calibration is recommended
Operation conditions	Temperature: 5°C-30°C, relative humidity up to 80%
Storage conditions	Temperature: 0°C-40°C, relative humidity up to 80%
Safety	EN 61010-1:2010+A1:2019, EN 61326-1:2013, EN IEC 61326-1:2021, EN 61000-6-1:2007, EN 61000-6-3:2007+A1
Warranty	24 months (calculated from the date of product activation or six months after the date of manufacture, whichever comes first)

Ordering Information

Name	Model	Order No.
BioAerosol Sampler-Economic I BAS-Eco	C2111	MACHC2111

Sampling Head	Order No. (Anodized aluminum)	Order No. (Stainless steel)
Φ 0.32mm*280 holes (28.3L/min)	J04A010140a	J04A010150a
Φ 0.60mm*280 holes (100L/min)	J04A010141a	J04A010151a
Φ 0.70mm*400 holes (200L/min)	J04A010142a	J04A010152a

MICRONVIEW

Tel: 1-520-528-1882

Email: mv@micronview.com

Website: <http://www.micronview.com>

Address: 2300 Highway 365 Suite 110, Nederland TX 77627, USA

©MICRONVIEW LIMITED

Version: 1.0

MICRONVIEW®

