


BAS

BioAerosol Sampler



Feature:

- Flow rate accuracy: $\pm 2.5\%$ (28.3L/min, 100L/min) / $\pm 3\%$ (200L/min)
- Long battery duration:
 - Battery last ≥ 10 hours @ 28.3L/min (16.98m³)
 - Battery last ≥ 10 hours @ 100L/min (60m³)
 - Battery last ≥ 4 hours @ 200L/min (48m³)
- Authority management (Divided into admin and operator)
- Dual-purpose sampler for general sampling and remote sampling (Supports up to 8m sample tubing)
- Supports various types of flowmeters for calibration, timesaving and low cost
- Complies with ISO 14698-1

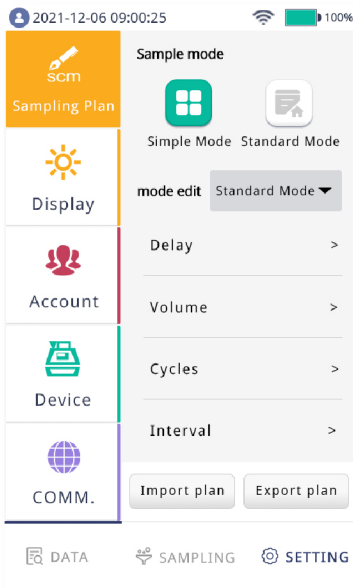

Setting



Smart and Easy

- Smart settings mode
- User-friendly input method

Data

Detailed and Flexible

- Intuitive sampling records
- Quick search capability for data records
- One-click for data exporting and printing

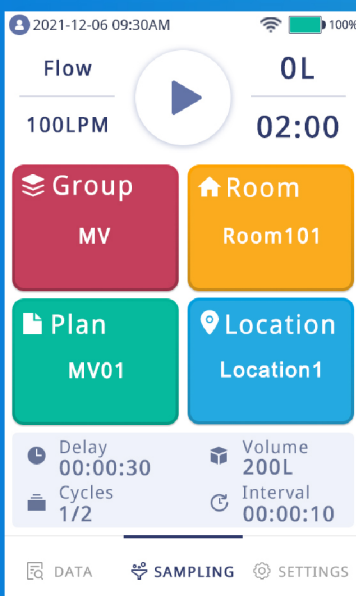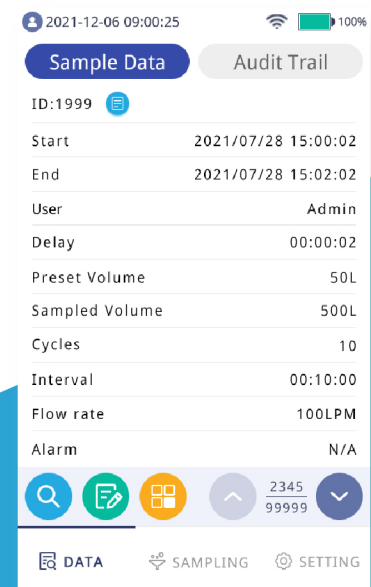

Sampling

Distinct and Direct


- Intelligent sampling mode enables continuous and cyclic sampling
- Clear and distinct sampling process
- LED lighting and beep prompt error

Isokinetic Sampling

When performing isokinetic sampling, the sampler is designed so that the airflow direction into the inlet of the sampling head is consistent with the average velocity of the unidirectional flow through the sampling head.



The flow velocity above the sampling head is 0.4m/s, consistent with isokinetic sampling.

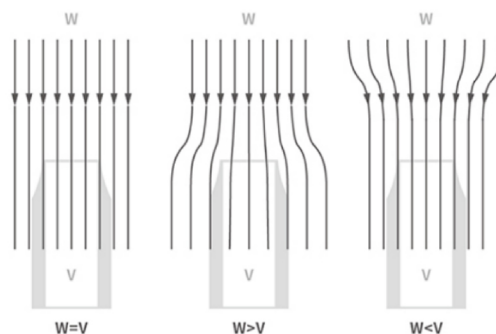


V=Sampling speed of sampling head
W=Speed of unidirectional flow

Remote Sampling

Used for sampling in an isolator, biological safety cabinet, horizontal laminar flow clean workbench and other significant zones.




Vertical flow sampling

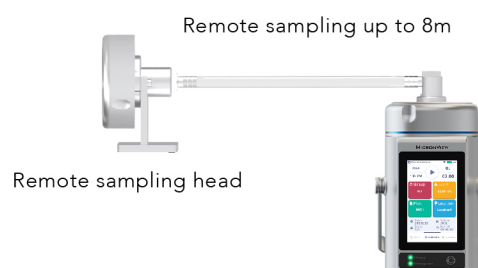

Horizontal flow sampling

Compressed Gas Sampling

Compressed gas is decompressed to atmospheric pressure by the High Pressure Diffuser (HPD) and transported to the BioAerosol Sampler for sampling.

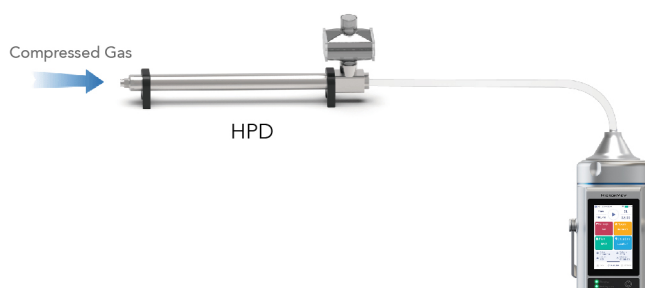

Specification Sheet

Specification	BioAerosol Sampler I BAS
Flow rate	28.3L/min $\pm 2.5\%$; 100L/min $\pm 2.5\%$; 200L/min $\pm 3\%$
Flow rate control	Built-in mass flow meter
Zone	Over 20000 locations
Sampling volumes	1-6000L, 6 default volumes
Delay	0-59 minutes 59 seconds
Cycles	1-99 cycles
Interval	0-59 minutes 59 seconds
Sampling head	$\Phi 0.32$ mm*280 holes (28.3L/min) $\Phi 0.60$ mm*280 holes (100L/min) $\Phi 0.85$ mm*280 holes (200L/min) (Anodized aluminum, 316L stainless steel material available, customizable sampling head)
Isokinetic probe caliber	$\Phi 101.8$ mm
Weight / Dimensions	2.7Kg (6.0lbs); $\Phi 110$ mm, H: 202mm (without foot mat)
External surface	Anodized aluminium
Power	AC 100-240V, 50 Hz/60 Hz
Battery	10.8V, 6800mAh rechargeable lithium battery Battery last ≥ 10 hours (28.3L/min) Battery last ≥ 10 hours (100L/min) Battery last ≥ 4 hours (200L/min)
Communication	USB, WIFI
Application extension	Connect with the Airborne Particle Counter (APC) and BioAerosol Monitoring System (BAMS) to collect contaminated samples in real time
Data storage	20,000 sets of records
Alarm	LED lighting and beep
Display	4.3" capacitive touch screen
Language	English, Chinese
Calibration frequency	Annual calibration is recommended
Operating conditions	Temperature: 5°C-30°C, relative humidity up to 80%
Storage conditions	Temperature: 0°C-40°C, relative humidity up to 80%
Safety	EN 61010-1:2010, EN 61326-1:2013, EN 61326-2-2:2013, EN 61000-6-1:2007, EN 61000-6-3:2007+A1
Warranty	2 years after activation

Ordering Information

Name	Model	Order No.
BioAerosol Sampler BAS	C130	MACHC130
BioAerosol Sampler (side handle) BAS	C230	MACHC230

Sampling Head	Order No. (Anodized aluminum)	Order No. (Stainless steel)
$\Phi 0.32$ mm*280 holes (28.3L/min)	J04A010140a	J04A010150a
$\Phi 0.60$ mm*280 holes (100L/min)	J04A010141a	J04A010151a
$\Phi 0.85$ mm*280 holes (200L/min)	J04A010142a	J04A010152a

MICRONVIEW

Tel: 1-520-528-1882

Email: mv@micronview.com

Website: <http://www.micronview.com>

Address: 2300 Highway 365 Suite 110, Nederland TX 77627, USA

©MICRONVIEW LIMITED

Version: 2.0

micronview.com

MICRONVIEW®

